

BEST AVAILABLE COPY




GB04/5310



INVESTOR IN PEOPLE

The Patent Office
Concept House
Cardiff Road
Newport
South Wales
NP10 8QQ

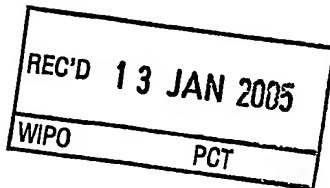

I, the undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.

In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in this certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L.C. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subjects the company to certain additional company law rules.



**PRIORITY
DOCUMENT**

SUBMITTED OR TRANSMITTED IN
COMPLIANCE WITH RULE 17.1(a) OR (b)

Signed

Stephen Hardley


Dated 8 November 2004



22 DEC 03 E960973 D02882
P01/7700-0.00-0329521.9 ACCOUNT/CHA

Request for grant of a patent

(See the notes on the back of this form. You can also get an explanatory leaflet from the Patent Office to help you fill in this form)

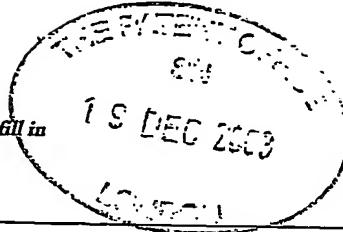

The Patent Office
Cardiff Road
Newport
South Wales
NP10 8QQ

1. Your reference

MJD/63952/000

0329521.9

2. Patent application number

(The Patent Office will fill this part in)

19 DEC 2003

3. Full name, address and postcode of the or of each applicant (*underline all surnames*)

KeyMed (Medical & Industrial Equipment) Ltd
KeyMed House
Stock Road
Southend-on-Sea
Essex SS2 5QH
United Kingdom

Patents ADP number (*if you know it*)

5486 93001

If the applicant is a corporate body, give the country/state of its incorporation

GB

4. Title of the invention A Dummy Medical Instrument for Use in a Simulator

5. Name of your agent (*if you have one*)

BOULT WADE TENNANT

"Address for service" in the United Kingdom to which all correspondence should be sent
(including the postcode)

VERULAM GARDENS
70 GRAY'S INN ROAD
LONDON
WC1X 8BT

Patents ADP number (*if you know it*)

42001

6. Priority: Complete this section if you are declaring priority from one or more earlier patent applications, filed in the last 12 months.

Country

Priority application number

(*if you know it*)

Date of filing

(*day / month / year*)

7. Divisionals, etc: Complete this section only if this application is a divisional application or resulted from an entitlement dispute (see note f)

Number of earlier UK application
(*day / month / year*)

8. Is a Patents Form 7/77 (Statement of inventorship and of right to grant of a patent) required in support of this request?

YES

Answer YES if:

- a) any applicant named in part 3 is not an inventor, or
- b) there is an inventor who is not named as an applicant, or
- c) any named applicant is a corporate body.

Otherwise answer NO (See note d)

Patents Form 1/77

9. Enter the number of sheets for any of the following items you are filing with this form.
Do not count copies of the same document

Continuation sheets of this form

Description 9

Claim(s) 2

Abstract

Drawing(s) 3 *X 3 SN*

10. If you are also filing any of the following, state how many against each item.

Priority documents

Translations of priority documents

Statement of inventorship and right to grant of a patent (*Patents Form 7/77*)

Request for preliminary examination and search (*Patents Form 9/77*) 1

Request for substantive examination (*Patents Form 10/77*) 1

Any other documents
(please specify)

11. I/We request the grant of a patent on the basis of this application.

Signature

Date

R. J. Wadsworth 19 December 2003

12. Name and daytime telephone number of person to contact in the United Kingdom Martyn J. Draper
020-7430-7500

Warning

After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked.

Notes

- a) If you need help to fill in this form or you have any questions, please contact the Patent Office on 08459 500505.
- b) Write your answers in capital letters using black ink or you may type them.
- c) If there is not enough space for all the relevant details on any part of this form, please continue on a separate sheet of paper and write "see continuation sheet" in the relevant part(s). Any continuation sheet should be attached to this form.
- d) If you have answered 'Yes' Patents Form 7/77 will need to be filed.
- e) Once you have filled in the form you must remember to sign and date it.
- f) For details of the fee and ways to pay please contact the Patent Office.

DUPPLICATE

- 1 -

A DUMMY MEDICAL INSTRUMENT FOR USE IN A SIMULATOR

The present invention relates to a dummy medical instrument for use in a simulator.

5

One type of simulator to which the present invention is applicable is that disclosed in GB A 2252656. This simulator simulates the operation of an endoscopic process.

10 A dummy endoscope is insertable into a fixture which is provided with a sensor mechanism to sense the longitudinal and rotational movement of the dummy endoscope. This information is fed to a controller which generates force feedback information based on virtual model data held in the computer memory. The force feedback applied to the dummy 15 endoscope is synchronised with a visual representation of the procedure so as to provide a realistic simulation providing a useful training tool to endoscope users.

With an instrument such as an endoscope, the tip of the 20 endoscope is manipulated by angulation control in the form of one or more control knobs on the handle of the endoscope which are linked to cables which extend down the insertion tube of the endoscope. Turning of the control knobs produces a corresponding movement of the cable and hence the 25 tip. An endoscope can have two control knobs one of which controls the left/right movement of the tip and the other of which controls the up/down movement of the tip.

An example of a dummy medical instrument which uses the 30 angulation cables to transmit the force feedback to the control body is disclosed in our earlier application WO 03/058583. In this instrument, the angulation cables are directed along the umbilical of the instrument to a main

unit where they are wound around motors which generate a variable force to provide force feedback.

In an operation such as a colonoscopy, the medical practitioner will often twist the control body a number of times to manipulate the endoscope along the colon. In doing this, the umbilical can become severely twisted. In order to remove the coils, the endoscopist or assistant may either reverse the control body at a suitable safe place in the colon, or unplug the umbilical and untwist it. This latter option is not possible with the arrangement disclosed in WO 03/658583 as the four wires extending from the umbilical are permanently wrapped around the force feedback motors in the main unit to which the umbilical is attached. Further, the permanent attachment of the umbilical to the main unit makes transportation and storage of the umbilical difficult.

To address this problem, we have proposed a dummy medical instrument with a connector which is detachable from the main unit by the user which contains means to apply a variable force to the cable. This is claimed in PCT/GB03/003106. This allows an operator to simulate the process which is carried out in practice of detaching and untwisting the umbilical before plugging it back into the main unit.

However, when the umbilical is detached in this way, sensing of the position of angulation controls is no longer possible. In practice, an operator will not generally move the angulation control when the umbilical is detached as his/her attention will be focused on untwisting the umbilical, rather than continuing to manipulate the endoscope. However, it is possible for the angulation

control to be manipulated at this time, and this will cause position to be lost.

According to the present invention there is provided a
5 dummy instrument for use in a simulator, the instrument comprising a control body with user manipulatable angulation control, an insertion tube and an umbilical extending from the control body, wherein in a real instrument corresponding to the one being simulated, at least one angulation cable
10 would extend from the user manipulatable controls to the tip of the insertion tube such that movement of the angulation control changes the angulation of the tip, and wherein in the dummy medical instrument the angulation cable extends from the user manipulatable angulation control, and down the
15 umbilical, the umbilical being releasably attached to a main unit, a motor within the instrument at the distal end of the umbilical to apply a variable force to the cable, and a position detector within the instrument to detect the angular position of the angulation control.

20

Thus, the umbilical is still detachable from the main unit. However, as the motor and position detector are within the instrument, when the umbilical is disconnected, the simulation will be aware of the position of the
25 angulation controls when the umbilical is reconnected. Also, as there is no detachable connection to the motors and sensors, there is no issue of having to re-engage with the motor and position detector when the umbilical is
30 reconnected, nor is there any danger of slippage across this releasable connection.

The position detector may either detect the position of the angulation control directly by measuring the rotation of the angulation control. Alternatively, it may detect this

indirectly measuring the position of the cable. In the latter case, the detector may be positioned at the distal end of the umbilical such that it can be packaged with the motor.

5

Preferably, the dummy instrument is provided with two pairs of angulation cables, each pair forming a loop around the control body, and around a respective motor at the distal end of the umbilical. Such an arrangement provides 10 the left/right and up/down control provided in a normal endoscope. Alternatively, the dummy instrument may be an instrument such as a bronchoscope which only has a single pair of angulation cables.

15

With use of the instrument, the cable loops can become slack. Therefore, preferably, means are provided to retension the loop.

As the end of the umbilical with one or more motors and 20 position sensors is intended to be removable, it should be designed such that the motor(s), in particular, are as light as possible. If the detached end of the umbilical is considered by users to be too cumbersome to manipulate, a connector at the distal end of the umbilical could be 25 configured to provide a two-part release, allowing release to a partially released position in which the umbilical may be rotated relative to the base unit, but in which the weight of the umbilical is still supported by the base unit, and a second fully released position in which the umbilical 30 is completely releasable from the base unit. Thus, if a user wishes only to untwist the umbilical, he can move the connector to its partially released configuration and untwist the umbilical without having to support the weight of the motors. If the instrument is to be replaced by a

different instrument, the connector is fully released and a replacement instrument can be introduced.

In order to provide full sensing of the angulation position, high resolution absolute position detecting is necessary. However, an absolute high resolution position detector suitable for this task is a large component. Therefore, preferably, the position of the or each angulation cable is sensed by a combination of a low resolution absolute position detector and a higher resolution incremental encoder. Such a combination provides the necessary absolute high resolution detection without requiring a single high resolution absolute position detector.

15

Examples of a dummy medical instrument constructed in accordance with the present invention will now be described with reference to the accompanying drawings, in which:

20 Fig. 1 is a plan view of the distal end of the umbilical and connector with the casing of the connector removed;

Fig. 2 is a perspective view of the arrangement of Fig. 25 1 from one end;

Fig. 3 is a perspective view of the arrangement of Fig. 30 1 from the opposite end.

The basic details of providing angulation force feedback by rerouting the angulation force feedback cables from the insertion portion of the endoscope back down the umbilical of the endoscope are described in detail in our

earlier applications WO 03/058583 and GB 0213981.4 and will not be described more fully here.

Suffice it to say that the four angulation cables in
5 the endoscope are routed around a set of pulleys in the instrument handle which turn them through 180° allowing them to be fed back down through the umbilical where they emerge at the distal end.

10 The four wires are grouped in pairs, namely a pair for controlling the up/down motion of the endoscope tip and a pair for controlling the left/right motion of the endoscope tip. Each pair is moved by rotation of one of the control knobs. In the dummy instrument, no actual movement of the 15 tip of the instrument occurs as the cables have been rerouted. However, the corresponding movement of the cables is sensed, and this information is fed into the simulation which detects the 'virtual' position of the endoscope tip and calculates the displayed view and the force feedback applied to the instrument accordingly.
20

The connector at the distal tip of the umbilical where the four cables emerge is shown in the accompanying drawings.

25 The distal end of the umbilical U is shown in the drawing and the four cables 1,2,3,4 are shown emerging from the distal end of the umbilical and passing into the connector.

30 The connector will have a cylindrical housing which has been removed from the drawings for clarity. This is connected between two circular end plates 5,6. End plates 5,6 are further supported by a pair of support plates 7,8

but Figs. 2 and 3 have been shown without the support plate 8 so as not to obscure the detail of the internal mechanism of the connector.

5 As best shown in Fig. 1, a number of spacers 9 external from the end of the umbilical U to the end plate 6. These spaces are covered by a moulded cover 10 which is not shown in Fig. 1. The purpose of these is to allow the space for the cables 1-4 to fan out to the extent shown in Fig. 3.

10

Each cable 1-4 consists of an inner wire 11 surrounded by a sheath 12. Each sheath 12 terminates at a bolt 13 which is screw threaded into a flange 14 in the housing. Such an arrangement allows the position of the sheath 12 to 15 be moved with respect to the wire 11, thereby allowing tensioning of the cable. Each wire 11 extends around a small pulley 15 which is mounted in a pulley block 16 and turns each wire through approximately 90° as best shown in Fig. 1. This angle may be different for each wire.

20

Also mounted within the connector are a pair of force feedback motors, namely an upper force feedback motor 20 and a lower force feedback motor 21. Each motor has an output shaft 22 to which a pair of pulley wheels 23,24 are rigidly 25 mounted.

In relation to the upper motor 20, the pulleys 23,24 are positioned such that they are aligned with the wires 11 from cables 1 and 2 respectively once these are routed 30 around the small pulleys 15. Specifically, the wire 11 from the first cable 1 is aligned with a groove in the top surface of the pulley 23, while the wire 11 from the second cable 2 is aligned with a groove in the lower surface of the

second pulley 24. The two wires are fastened to the two pulleys by respective screw clamps.

A similar arrangement is employed in respect of the
5 lower motor 21 and the third and fourth cables 3,4. Thus,
the force output from the two motors can be transmitted to
the pairs of cables 1,2 and 3,4 to provide force feedback
for the up/down and left/right movements of the endoscope
tip respectively.

10 Position readings for the cable wires are provided by a
combination of an incremental (relative) position detector
30 and an absolute position detector on board 31 for each
motor. An absolute position detector is required as the
15 connection is not powered when disconnected from its unit.
However, a high resolution absolute position detector would
be too large to fit in the space available. A combination
of the absolute position detector with a high resolution
incremental encoder 30 overcomes this problem.

20 A number of features are provided on the end plate 5 to
interface with the base unit when the connector is in place.
In particular, an alignment orifice 40 is arranged to
cooperate with a corresponding pin on the base unit to
25 ensure that the connector is inserted in the correct
orientation. Two air ports 41 are provided. These connect
to an air supply in the base unit, and tubes run from these
ports along the connector and up the umbilical to the
endoscope handle. One of these provides the air to the
30 insufflation button on the handle which provides a realistic
feel to this button, while the other provides a supply of
air for a pneumatic force feedback device. A power input
port 42 is arranged to provide a connection for power to the

motors 20,21 and a data board 42 provides a connection for data to and from the connector.

CLAIMS

1. A dummy instrument for use in a simulator, the instrument comprising a control body with user manipulatable angulation control, an insertion tube and an umbilical extending from the control body, wherein in a real instrument corresponding to the one being simulated, at least one angulation cable would extend from the user manipulatable controls to the tip of the insertion tube such that movement of the angulation control changes the angulation of the tip, and wherein in the dummy medical instrument the angulation cable extends from the user manipulatable angulation control, and down the umbilical, the umbilical being releasably attached to a main unit, a motor within the instrument at the distal end of the umbilical to apply a variable force to the cable, and a position detector within the instrument to detect the angular position of the angulation control.

2. An instrument according to claim 1, wherein the position detector measures the rotation of the control.

3. An instrument according to claim 1, wherein the position detector measures the displacement of the cable.

4. An instrument according to claim 3, wherein the position detector is located at the distal end of the umbilical.

5. An instrument according to any preceding claim, wherein the instrument is provided with two pairs of angulation cables each pair forming a loop around the control body, and around a respective motor at the distal end of the umbilical.

6. An instrument according to claim 5, wherein means are provided to retension the loop.


5 7. An instrument according to any one of the preceding claims, wherein a connector at the distal end of the umbilical is configured to provide a two-part release, allowing release to a partially released position in which the umbilical may be rotated relative to the base unit, but
10 in which the weight of the umbilical is still supported by the base unit, and a second fully released position in which the umbilical is completely releasable from the base unit.

15 8. An instrument according to any one of the preceding claims, wherein the position of the or each angulation cable is sensed by a combination of a low resolution absolute position detector and a higher resolution incremental encoder.

20

25

#492671 v5

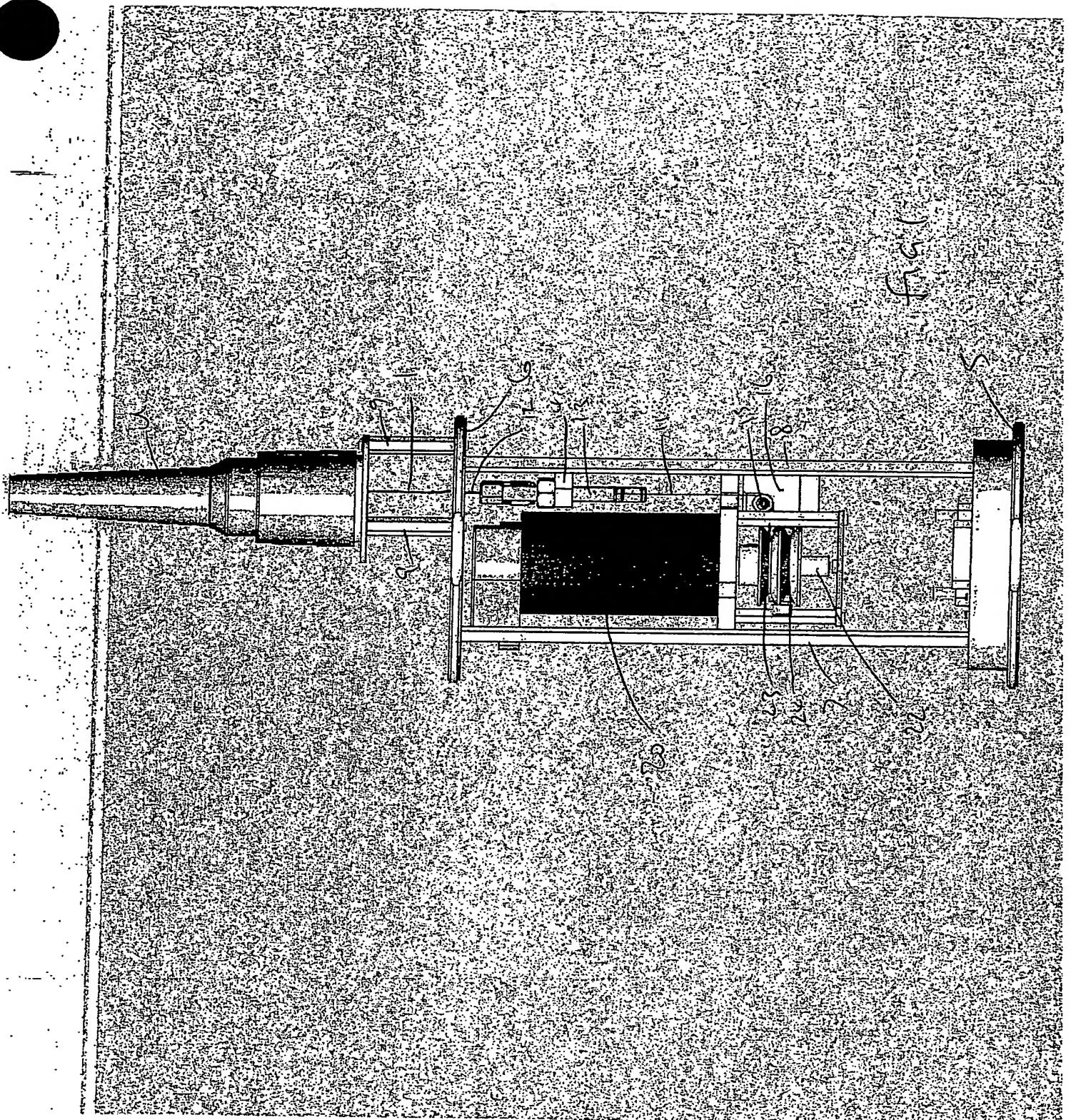


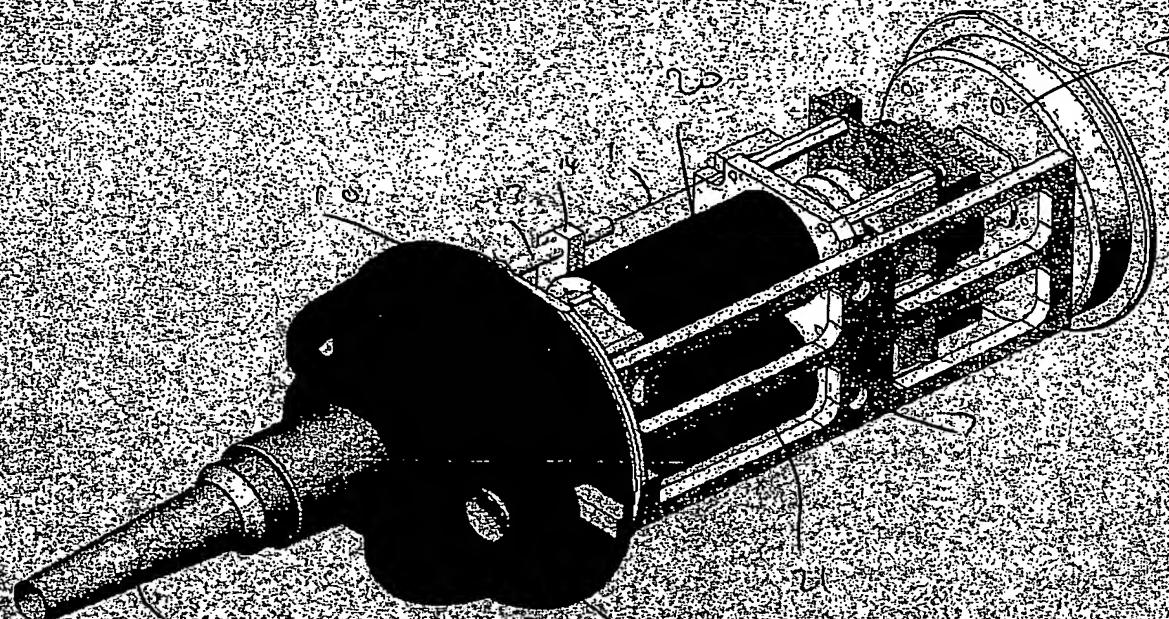


Fig 2

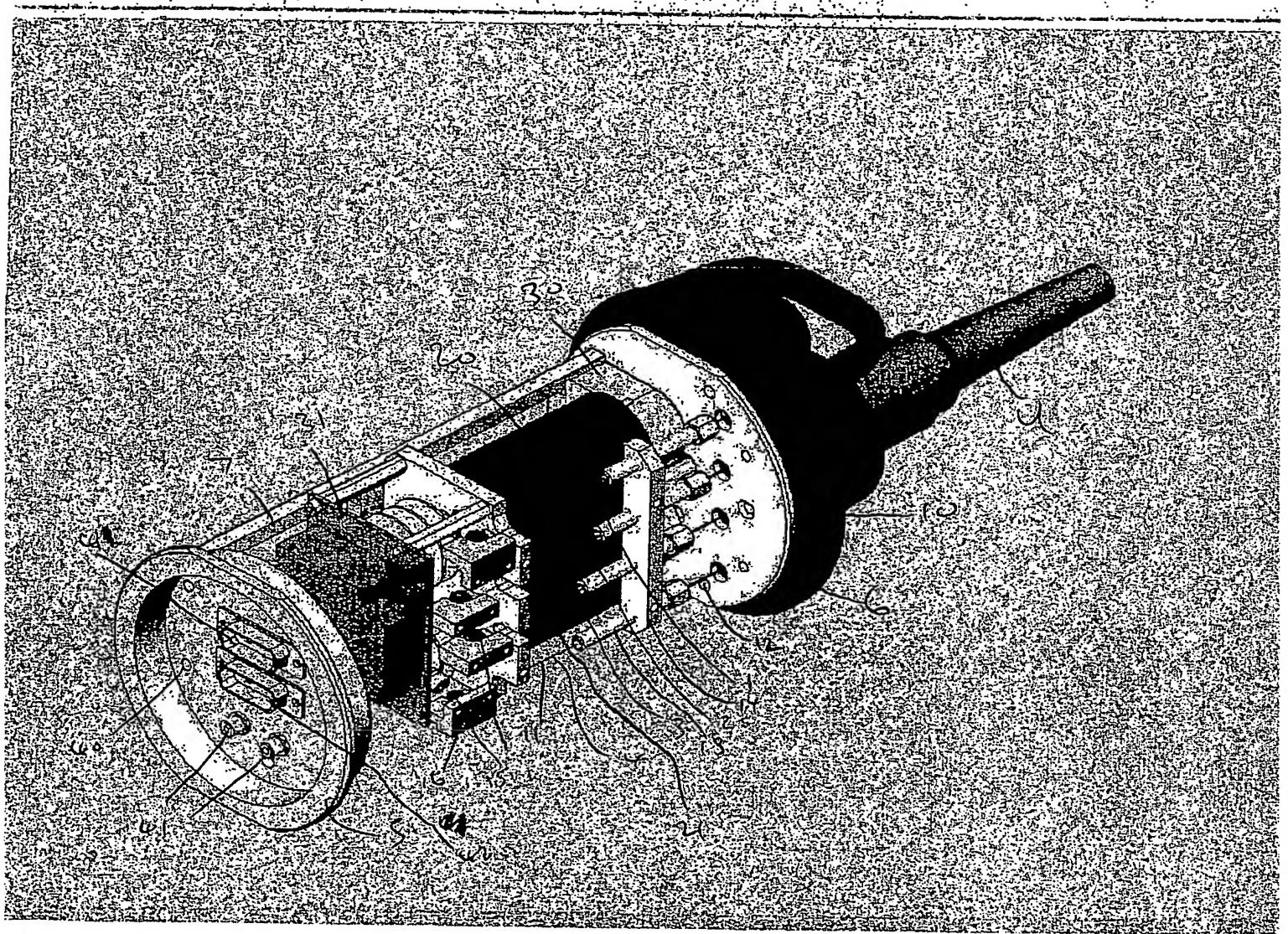

F C 3

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record.**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- BLACK BORDERS**
- IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- FADED TEXT OR DRAWING**
- BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- SKEWED/SLANTED IMAGES**
- COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- GRAY SCALE DOCUMENTS**
- LINES OR MARKS ON ORIGINAL DOCUMENT**
- REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- OTHER: _____**

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.